


## Jain Laboratory Instruments





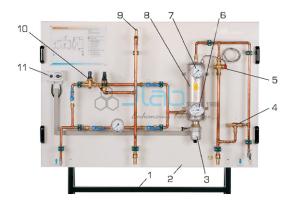

**Product Code. JL-ELA-5449** 

## **Safety Devices Training Panel**

## **Description**

Safety Devices Training Panel

An electrically heated pressure vessel is mounted on the panel.

The vessel can be filled via different pipe sections.

One pipe section is equipped with a system of safety devices for protecting a sealed water heater.

The pressure vessel is also equipped with a thermal discharge safety device that discharges to a blow-off pipe.

When the pressure or temperature level set is exceeded, the safety devices are tripped.

Pressure and temperature in the system are displayed.

A supply of water is required for operation.

The training panel can be fitted to the Universal Stand.

Features:

Trainer on heating systems and plumbing

Function of safety devices

Tested pressure vessel

